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CONTENTS

CHAPTER 1 START

1.1 DKMC01 PACKAGE INCLUDES	2
1.2 THREE WAYS TO INSTALL DKMC01 TO MICROSCOPE	2
1.3 ANYTY MICROSCOPE SOFTWARE INSTALLATION INSTRUCTIONS	3

CHAPTER 2 INTRODUCTION OF ANYTY MICROSCOPE SOFTWARE

2.1 ANYTY MICROSCOPE SOFTWARE PANEL	6
2.2 OPERATE LANGUAGES	8

CHAPTER 3 PICTURE PREVIEW

3.1 PICTURE PREVIEW WINDOW	9
3.2 PICTURE CONTROL BAR	10
3.3 PICTURE PARAMETERS CONTROL BAR	10

CHAPTER 4 SNAPSHOT AND RECORD VIDEO

4.1 SNAPSHOT	11
4.2 RECORD VIDEO	11

CHAPTER 5 PICTURE AND FILE MANAGEMENT

5.1 PICTURE AND VIDEO FILES BROWSE BAR	12
5.2 PICTURES BROWSE DIRECTORY	13
5.3 PICTURE FILES WINDOW	13
5.4 EDIT AND SAVE PICTURE/FILE, PDF CONVERSE	14
5.5 DISPLAY VIDEO	15

CHAPTER 6 : WARRANTY CARD AND OTHERS

6.1 WARRANTY CARD	16
6.2 ANYTY MICROSCOPE SOFTWARE UPDATE SUPPORT	17

Chapter 1 Start

1.1 DKMC01 Package Includes:

- A、DKMC01 body
- B、0.5 × minifier (23.2mm)
- C、30 mm adapter
- D、30.5 mm adapter
- E、operate software CD
- F、certification

1.2 Three Ways to Install DKMC01 to Microscope

Please choose suitable way according to microscope interface.

First install way:

Second install way:


Third install way:

1.3 Anyty Microscope Software Installation Instructions

Step 1: Start the computer and put the installation disk into the CD-ROM.


Step 2: Select " Anyty Microscope Install.exe "file. When the following window appears.

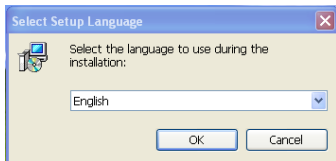
Please choose the language and click 

Tip: Please ensure you have installed correspondent language pack on your computer; otherwise, messy code would appear.

Step 3: When the window "Welcome to the Anyty Microscope Setup Wizard" appears, please click. 

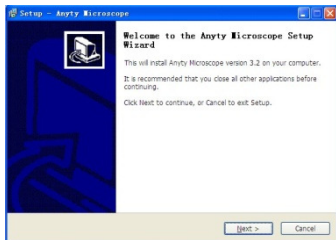Step 4: When the following window "Select Destination Location" appears, click "Browse" to choose the file installation path. Usually, you can keep clicking "next" if you don't need to change any information for installation.

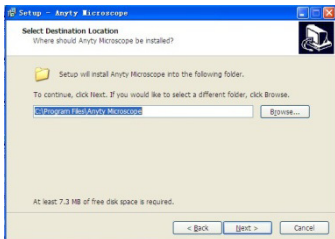Step 5: When the following window appears, the software is installed successfully.


Click  to complete the installation.

Step 6: Plug the 3R-DKMC01's USB cable into computer. The "Found New Hardware"

dialog box appears, computer system will automatically complete installation of driver program.

Step7. After the window closes, check the Anyty Microscope icon and start it on your desktop. If you cannot find it, check "Start" → "All Programs"→" Anyty Microscope ".

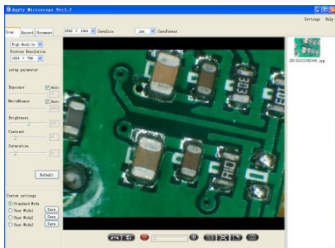
Tip 1: Before starting "Anyty Microscope", please make sure 3R-DKMC01 has connected to computer properly. If the device is connect to more than one USB port for the first time, the system might prompt you to install driver program, the system will automatically complete the driver installation for every USB port.

Tip 2: Please close " Anyty Microscope " program first before unplug 3R-DKMC01's USB cable from the device.

Chapter 2 Introduction of Anyty Microscope Software

2.1 Anyty Microscope Software Panel

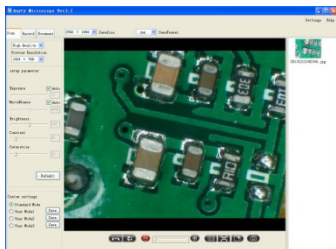Resolution Settings

Picture resolution: Anyty Microscope software allows users to choose different preview and picture resolution listed below :

640×480,800×600,1280×720,1024×768 and 2592×1944


Preview resolution influent only preview image quality and speed but not picture quality.

If you choose 800×600, Anyty Microscope software will take picture of image on preview window; if you choose 2592×1944, Anyty Microscope software will take picture of what it can fully captured.


If you have zoom in/zoom out the image, taken pictures might be different from preview image. Default resolution is 2592×1944.

Video resolution: Picture resolution has 640×480 , 800×600,1024×768,1280×720,1280×960,and 2592×1944 for selection.

The higher resolution, the better effect of video, and the larger storage space will need. Default video resolution is 1024×768.

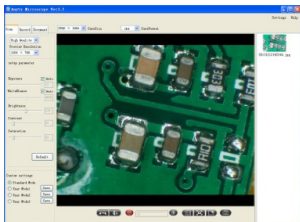Picture save path: click “Settings” on top right panel corner to show below manual, choose save path according to your need.

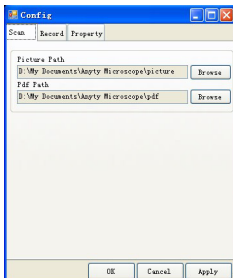
Tip 1: Default save path is (D:\My Documents\Anyty Microscope\picture)

Tip 2: Please add a new folder in computer operate system beforehand if you need, it don't support adding new folder on pop-up browse window.

Tip 3: Default save format is JPG, please choose other format if you need.

Video save path: click "Settings" on top right corner to show below manual, choose save path according to your need.

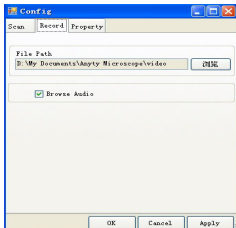
Tip 1: Default save path is (D:\My Documents\Anyty Microscope\video)

Tip 2: Please add a new folder in operate system beforehand if you need, it don't support adding new folder on pop-up browse window.


2.2 Operate Languages

You may choose operate language by clicking " settings" on the right top of the software panel: English, simplified Chinese, Traditional Chinese and Japanese.

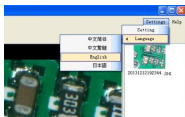Chapter 3 Picture Preview

Picture preview operate interface includes three parts: picture preview window, picture control bar, and picture parameters control bar.

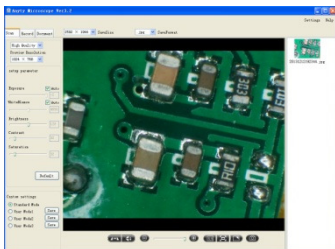3.1 Picture Preview Window

Real-time display observe picture, this window used for checking observed picture.

In order to get better display effect, you can click the full screen icon on the bottom-right of the software panel.

Right-click your mouse or press Esc key to exit full screen display


3.2 Picture Control Bar

Zoom in/zoom out: click "+" and "-" button, or click the seek bar directly with your mouse to zoom in and out the preview screen. It has six levels of zoom in/out. Default to intermediate position with smallest preview screen.

3.3Picture Parameters Control Bar

In order to obtain better display effect, you can adjust exposure value, brightness value, contrast value, and degree of saturation according to observed object and personal preference. It will back to default settings when you restart the software. If you need to save revised parameters, please refer to method below.

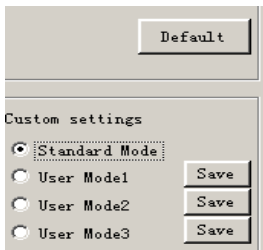Save parameter settings: you can revise and save the parameters according to personal preference. Methods: Click "save" in "custom define" after revising the parameters.

Chapter 4 Snapshot and Record Video

4.1 Snapshot

Click  button to snapshot

Tip 1: snapshot resolution will be determined by picture resolution settings.

Tip 2: picture save path default to software installation path, please click "save as" if you need other save path.

4.2 Record Video

Click  button to start to record video

Click  button to stop recording video

Tip1: video resolution will be determined by video resolution settings.

Tip2: video save path default to software installation path, please click "Settings"→"Record"→"file path" on the top right corner of the software panel if you need other save path.

Tip 3: it will need more storage space for recording, please ensure your have enough hard disk space for it.


Chapter 5 Picture and File Management

5.1 Picture and Video Files Browse Bar

Open picture files: browse pictures on picture folder with your computer's file explore;

Picture saves as: save picture with selected path, name, and format.

Picture formats: pdf, bmp, jpg, png.

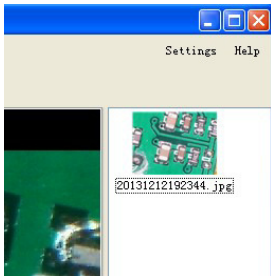Copy and paste picture to paste board: you may copy and paste picture to computer's paste board for other programs.

Pictures zoom in and zoom out: zoom in and zoom out the preview picture; click "actual size" will display picture with original resolution; click "adapt screen" will display picture with resolution suitable for preview window.

Back to default picture save directory: you can easily back to default directory by clicking this button if you have chosen other directory to open pictures before.

Back to default video save directory: you can easily back to default directory by clicking this button if you have chosen other directory to display video before.

5.2 Pictures Browse Directory

D:\MyDocuments\Anyty Microscope\picture

Used for choosing picture browse directory.


5.3 Picture Files Window

Here shows all picture thumbnails on browse directory, double-click thumbnail to display on preview window. You may rename pictures by right-clicking it under picture selected status.

5.4 Edit and Save Picture/File, PDF Converse

Click “Doc”, you may crop, screen shot and rotate the picture. Show as below.

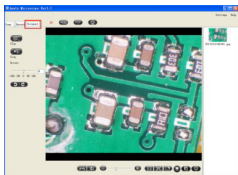
Export Picture

Click “Export”, select save path, file name, and click “save”

PDF format converse: click “PDF” to choose save path, input file name, and click “save”.

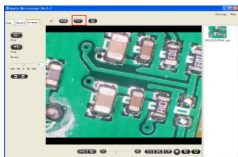5.5 Display Video

Open or double-click “back to video fault directory” button to pen video save directory, double-click the video file to display it with computer default video player.

Tip: please make sure your computer operate system has installed video player support AVI format.

Chapter 6: Warranty Card and Others

6.1 Warranty Card

Product	Anyty Microscope		
model	3R-DKMC01		
Serial NO.			
Sale Date			
User name:		Mobile/Tel	
Distributor's Tel			
Stamp of distributor			
Valid only with stamp of distributor			

To ensure your interests, please ask distributor to write down sales store name, address, date of purchase and stamp on your warranty card.

Warranty valid date is one year from the date of purchase

Guaranteed service: free repair for one year if hardware failure caused by producing.

Non-guaranteed service: conditions mentioned below are not including on warranty service:

- 1, Use improperly, detach or overhaul or modify by yourself (or **commit to other**)
 2. Failure or damage due to acts of the nature, for example, nature disaster, lightning, abnormal electronic voltage, and the occurrence of environmental factors
 3. Serial Number wrong, warranty card content fake, wrong, or too blur to identify.
- Service below shall be charged when warranty date expires: 1. Product test fee (includes shipping cost) 2. Maintenance fee 3. New component fee.

6.2 Anyty Microscope software Update Support

3R group will update the software of Anyty series microscope products occasionally, please download the latest version from the products page on our company's website if you need.

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